




Pre-shaded Multi with Smart Incisal™

- Multi Chroma and Blue Incisal Enhancer Features Incorporated Into One Disc
- Anterior Translucency* (48%, STML)
- Maximum Strength (>1200 MPa**)

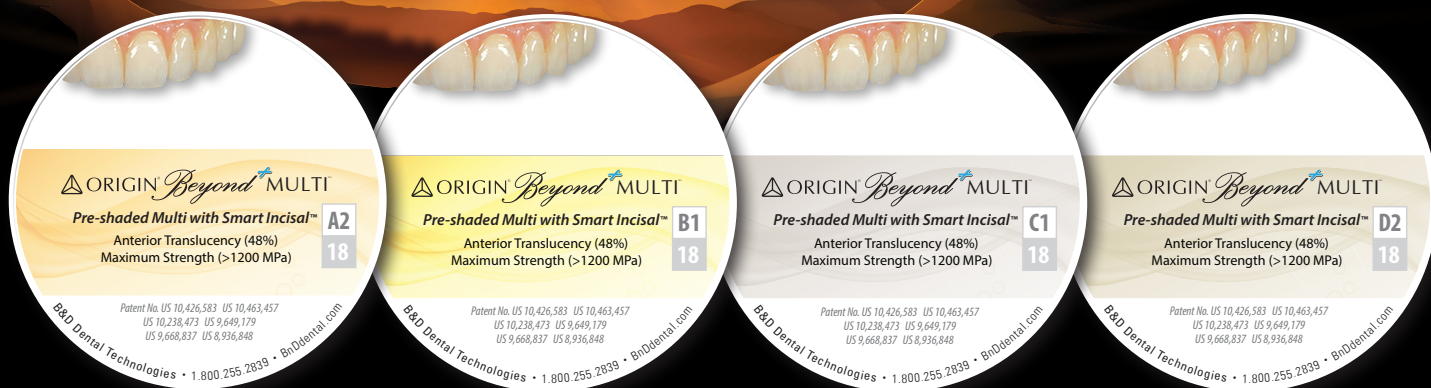

High Performance
Full Strength/Super Translucency
Ideal for Both Anterior & Posterior



Beyond+ Multi™ - layered discs have **MORE INTENSE CHROMA IN THE CERVICAL** and **DECREASING CHROMA TOWARDS THE INCISAL AREA**, offering the **MOST EFFICIENT PRODUCTION** speed combined with the natural look found in the shade tab. Each disc will give you the **EXACT SHADE** after simply milling and sintering. In addition, the incisal blue gray effect has been incorporated on the top 1/3 area of the disc to **REDUCE THE WHITE L* VALUE** and **PRESENT THE NATURAL TOOTH BLUEISH GRAYISH TRANSLUCENCY**.

(Patented Technology - Patent No. US 10,426,583 US 10,463,457 US 10,238,473 US 9,649,179 US 9,668,837 US 8,936,848)

*Total Light Transmittance level: minimum 48% at 600 nm wavelength, PerkinElmer Lambda 35 UV/VIS Spectrometer, Evaluated with raw material (white color).

**3 point bending strength according to ISO 6872:2015. Evaluated with raw material (white color).

**1100 MPa : up to 14 mm, 1200 MPa: over 16 mm

1. Anterior Translucency | Posterior Strength


- ◆ Higher translucency in the incisal area and higher chroma in the cervical area
- ◆ Decreasing chroma towards the incisal area
- ◆ Efficient production speed
- ◆ A natural gradation comparable to the shade tab

2. Smart Incisal™

EXACT SHADES – SIMPLY MILL & SINTER


Patented **Smart Incisal™** technology incorporated in the top 1/3 area of the disc.

- ◆ Mimics natural and subtle blueish gray translucency
- ◆ Increased incisal effect

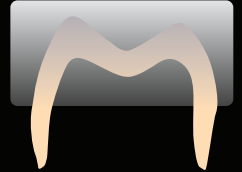


3. Unrivalled Strength

Multilayer Zirconia for Monolithic Restoration for Both *Anterior & Posterior*

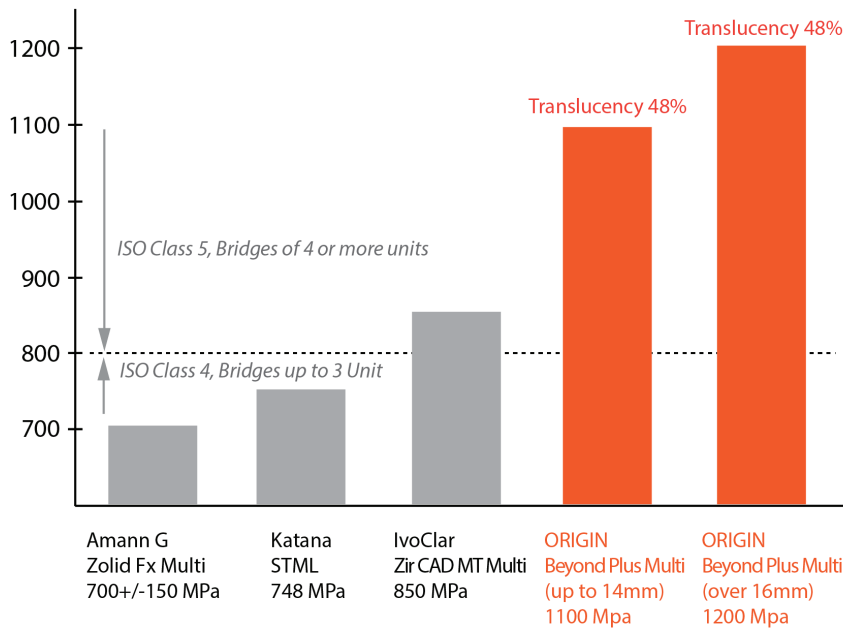

ORIGIN Uses Different Technology - Colloidal Pressing

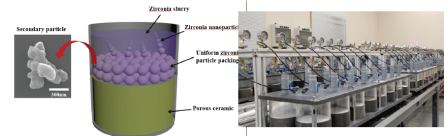

Secondary Particle Size: 300 nm



Industry's lowest shrinkage rate

AVAILABLE IN ALL 16 CLASSIC SHADES + BLEACH SHADE

For those who prefer to keep a lower inventory, customized liquids provide easy conversion to less common shades.

ex. Easily convert A2 disc to A3, A3.5 or A4 crowns.

CHROMA™
Version 8.1

